


ZoneOne Stability Test Chambers

Light-duty temperature/humidity chambers featuring a large volume and full-view doors

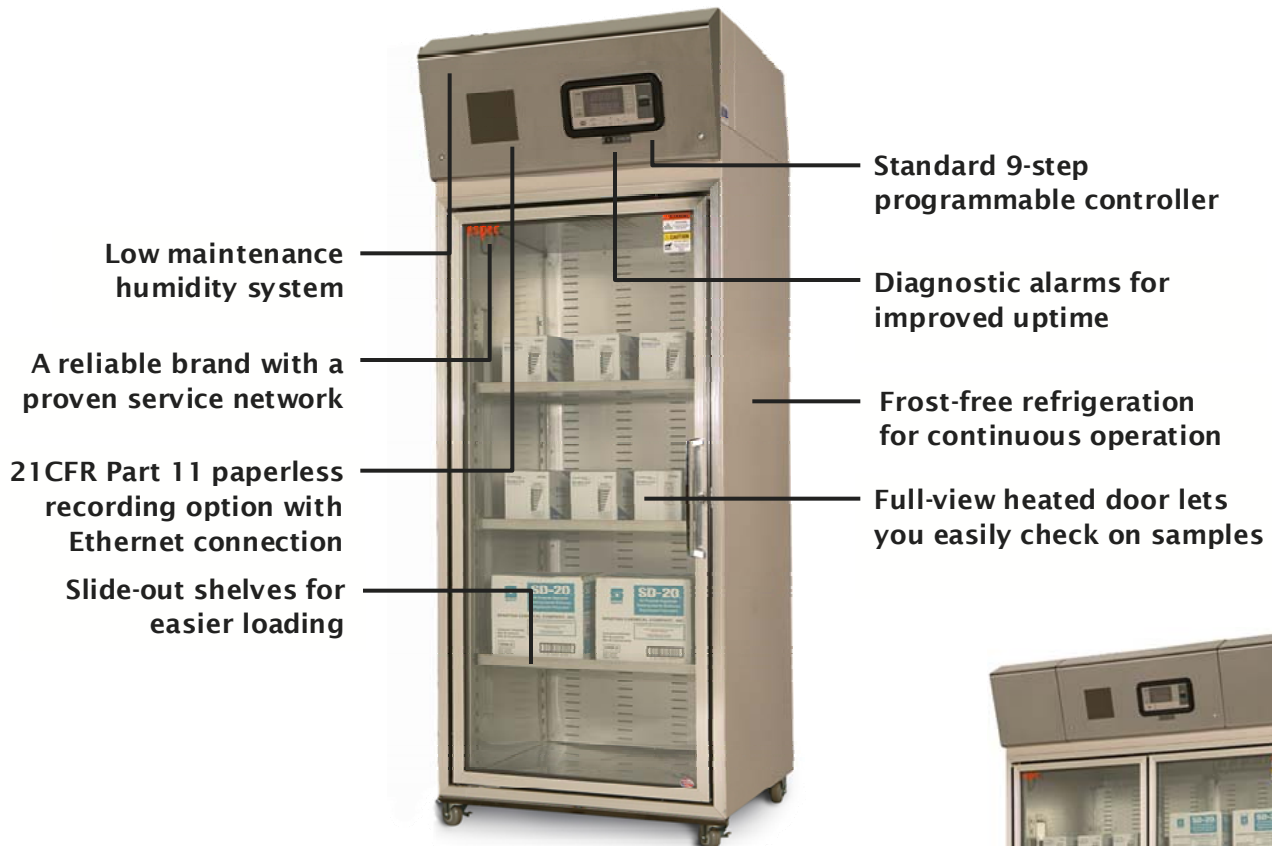


New double-door model

Features

- Performance meets ICH guidelines for pharmaceutical stability testing
- Two sizes available: 28 or 56 cubic foot volume
- Heated full-view window for easy viewing without condensation problems
- Horizontal airflow for uniform sample exposure
- Easy-to-use controller with nine-step programming for automatic start/stop and ramping
- Includes three adjustable solid shelves with lips to protect against spills
- Designed for flush installation on all sides; casters are standard
- One year parts and labor domestic warranty with factory-trained technicians
- IQ OQ PQ certification services available

ZoneOne Specifications

Model	EZO30-EA (single door)	EZO60-EA (double door)
Interior Volume	28 cu. ft. (800 L)	56 cu. ft. (1,600 L)
Interior Dimensions	27.5" x 27.8" x 63.5" (70 x 71 x 161 cm)	55" x 27.8" x 63.5" (140 x 71 x 161 cm)
Exterior Dimensions	34" x 40" x 91.5" (87 x 102 x 233 cm)	67" x 40" x 91.5" (172 x 102 x 233 cm)
Temperature Range	10 to 65°C (50 to 150°F)	
Temperature Fluctuation	±0.25°C (between 20 and 40°C)	
Humidity Range	20 to 90% RH (see chart)	
Humidity Fluctuation	±2.0% RH, above 20°C	
Power Required	208/230V, single phase, 60Hz	
Water Required	De-ionized water ¼" compression fitting	
Gravity Drain	½" hose fitting	

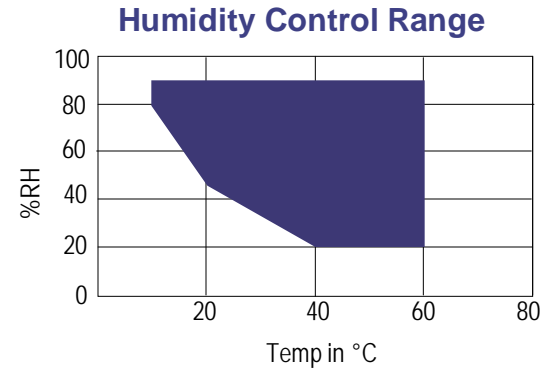

Programmer/Controller

- Designed specifically for temperature/humidity test chambers
- Constant mode for stability testing
- Nine-step program mode for automated ramping and start/stop
- Large LED display
- Alarms for out-of-setting conditions (temperature and humidity)
- Redundant overheat protection
- Includes either RS-485, RS-232C, or GP-IB computer interface
- LabView drivers for GP-IB or ER-100SII software for RS-232



Options Available

- Strip chart or circular chart recorders
- Additional shelves (three are included)
- Water filter cartridge system
- Water supply tank
- Cable ports for running-in wires: 25mm or 50mm diameter
- Interior light
- Stainless steel exterior
- Power cord

Option: Paperless Recorder

- ◆ Secure data storage complies with CFR21 Part 11 requirements
- ◆ Color LCD display
- ◆ Records to CompactFlash
- ◆ Ethernet connection for intranet access and datalogging
- ◆ Alarms can trigger email alerts

ESPEC NORTH AMERICA
4141 Central Parkway
Hudsonville, MI 49426
Phone: 616-896-6100
www.espec.com